



HYPERLITE 1 PORTABLE HYPERBARIC CHAMBER



PORTABLE | FOLDABLE | INDISPENSABLE

- ✓ The original flexible hyperbaric chamber
  - ✓ Internationally approved
  - ✓ Lightweight & easily transportable
  - ✓ Performs all 100% oxygen therapies
  - ✓ Multiple lifesaving applications
  - ✓ In service in over 30 countries

*“Now the ability exists to rapidly deliver acute and effective hyperbaric treatment even in remote locations”*

International Hyperbaric Medical Foundation




Delivering on-scene hyperbaric oxygen without delay, the Hyperlite 1 is the only truly portable non-metallic system capable of all 100% oxygen treatments either on-scene or during transfer by land, air or sea.

# HYPERLITE 1

## PORTABLE HYPERBARIC CHAMBER

Over 150 units in service in over 30 countries, the Hyperlite 1 was the first commercially available non-metallic chamber and remains the industry leader for military, commercial and sport diving industries worldwide.

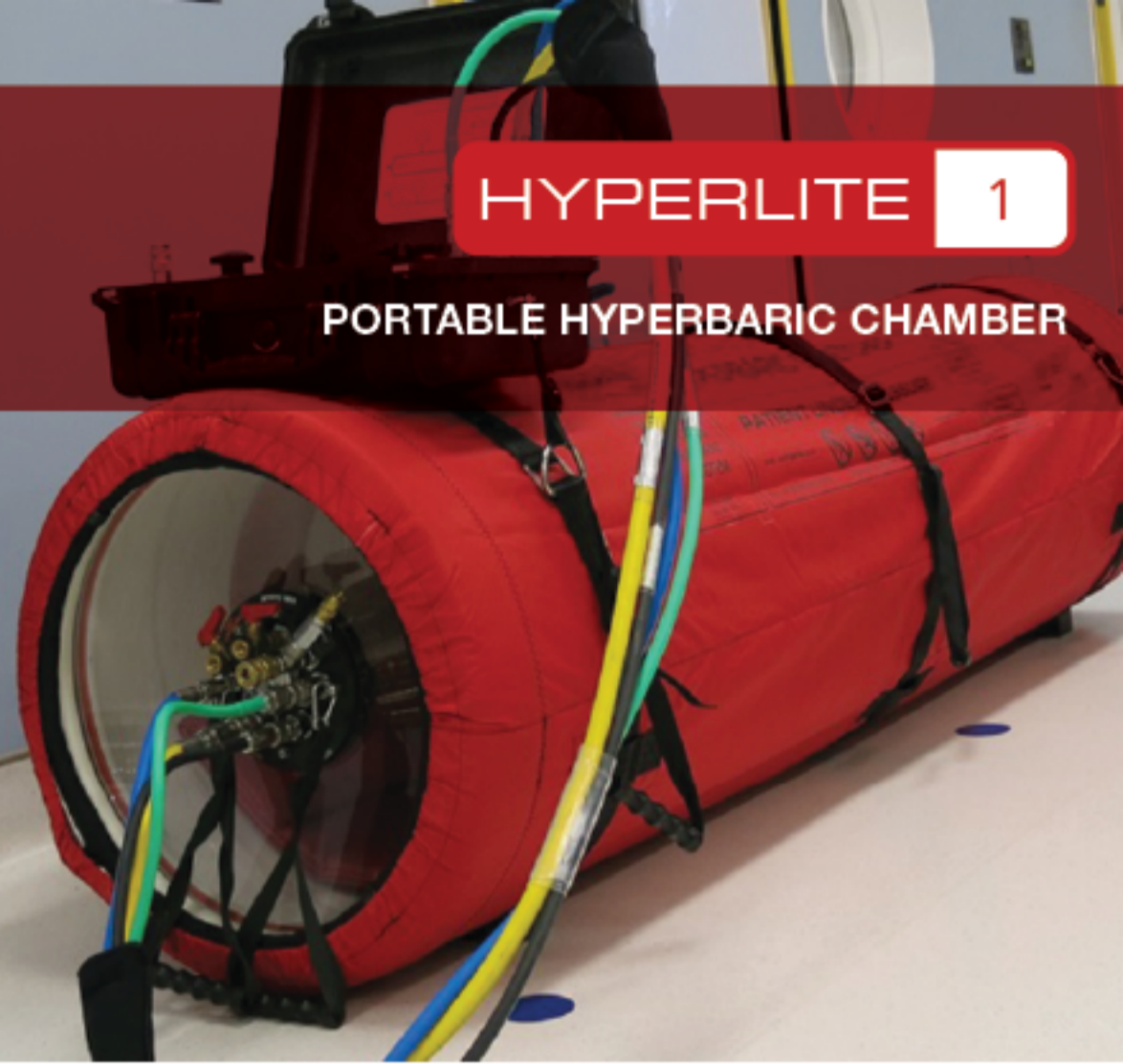

### KEY FEATURES

Pressure Vessel Weight:	50kg (110 lb)	Operating Pressure:	2.3 Bar(g)/23msw (75fsw/33.5psi(g))
Weight in box:	100kg (229 lb)	Operational Length:	224.5cm (88.5")
Box Dimensions (3):	68cm x 68cm x 64cm (27" x 27" x 25")	Tube Diameter:	59.5cm (23.5")
	47cm x 36cm x 18cm (19" x 14" x 7")	System life:	Up to 30 years (system upgrade available every 10 years)
	47cm x 36cm x 18cm (19" x 14" x 7")		



A vital first response solution for injured divers, high altitude pilots, submariners and carbon monoxide sufferers. Treatment within the first hour dramatically reduces the risk of serious injury or even death.

### SYSTEM BENEFITS



- Capable of all 100% oxygen treatment tables
- Continued treatment during transportation
- Small & lightweight packaged volume – Where space is at a premium
- Treatment up to 5,500 (18,000 ft) in an unpressurised cabin
- Minimal gas requirements – Uses 15 litre SCUBA Cylinders

Worldwide Approvals including: -  
 ASME/PVHO-1 | EC Approval | FDA Complaint (510k) | ISO: 13485  
 Approved Quality Management System | Health Canada Approved |  
 Lloyds Register witnessed testing | US Navy Diving & Manned  
 Hyperbaric Systems Safety Certification

- Suitable for helicopter winch
- Simple control system



For more information on the Hyperlite 1 system, our comprehensive operational & medical training courses, finance options and other first response medical solutions, contact us today at [www.sosgroup.co](http://www.sosgroup.co)